

**COUNTERTOPS:**

The surfaces below are for fume hood work areas along with standard countertop space within the lab.

The countertops can be ordered in various lengths and sizes with added features depending on your needs. Sinks along with fixtures can be supplied with the countertops. The countertops come in two forms epoxy resin which can be dished or flat depending on your needs. The resin countertops are 1 1/4" thick and when there is a dish it is 1/4" deep. The epoxy resin countertops are a special order item with standard sizes and sink cutout locations. The diagram on the next page shows the standard cup sink locations for a 48 CA style hood.

We also offer polished 304 stainless steel counters dished or flat along with the option for sinks and back splashes. The stainless steel tops are constructed of no less than 16 gauge material with internal reinforcing. All areas that are exposed should have a polished finish. The stainless steel countertop is built to your needed lengths and depths. The standard height is 3/4" to 1 1/4" and when there is a dish it is 1/4" deep typically. Also standard sink sizes are ordered as needed based on the customers needs.


Typical cup-sink size will 3" x 6" with a strainer located in the areas described on the next page or placed according to specifications. The sink is fabricated from plastic unless specified to be of another material, color is black. Other sink sizes are available upon request.

See the chart below for more information:

<b>COUNTERTOP SIZES FOR FUME HOODS</b>				
STOCK COUNTERTOPS ARE ALL DISHED 1/4" DEEP				
PART NO.	LENGTH	DEPTH	THICKNESS	MATERIAL
905-00S-035	35"	30"	1 1/4"	304 STN. STL.
905-00S-036	36"	30"	1 1/4"	304 STN. STL.
905-00S-047	47"	30"	1 1/4"	304 STN. STL.
905-00S-048	48"	30"	1 1/4"	304 STN. STL.
905-00S-058	58"	30"	1 1/4"	304 STN. STL.
905-00S-060	60"	30"	1 1/4"	304 STN. STL.
905-00S-070	70"	30"	1 1/4"	304 STN. STL.
905-00S-072	72"	30"	1 1/4"	304 STN. STL.
905-00S-094	94"	30"	1 1/4"	304 STN. STL.
905-00S-096	96"	30"	1 1/4"	304 STN. STL.
905-00R-035	35"	30"	1 1/4"	EPOXY RESIN
905-00R-036	36"	30"	1 1/4"	EPOXY RESIN
905-00R-047	47"	30"	1 1/4"	EPOXY RESIN
905-00R-048	48"	30"	1 1/4"	EPOXY RESIN
905-00R-058	58"	30"	1 1/4"	EPOXY RESIN
905-00R-060	60"	30"	1 1/4"	EPOXY RESIN
905-00R-070	70"	30"	1 1/4"	EPOXY RESIN
905-00R-072	72"	30"	1 1/4"	EPOXY RESIN
905-00R-094	94"	30"	1 1/4"	EPOXY RESIN
905-00R-096	96"	30"	1 1/4"	EPOXY RESIN

**STAINLESS STEEL –**

Counter tops are fabricated from type 304 stainless steel a minimum of 16 gauge. All edges are ground buffed and polished to a number 4 finish. All tops use hat channels and added support with sound absorbing material for quiet use. Contact your dealer for more details, will quote all shapes and sizes to help fill your requirements.

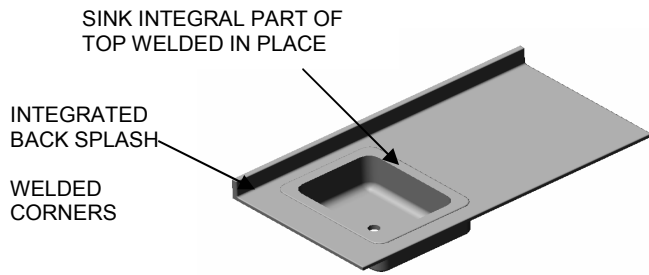

**A STAINLESS STEEL COUNTERTOP WITH A SINK AND BACK SPLASH.**

**COUNTERTOPS**

